SINCE 1908

*Solis*

# THERMOPAD

Typ/Type/Tipo 2260

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
## ▲ DESCRIPTION OF APPLIANCE PARTS

- A Heating pad, consists of heating element and removable protective cover
- B Control unit with sliding switch for switching the appliance on and off and for selecting one of 3 temperature levels; with temperature level display
- C Power cord with power plug and connection adapter

**GB**

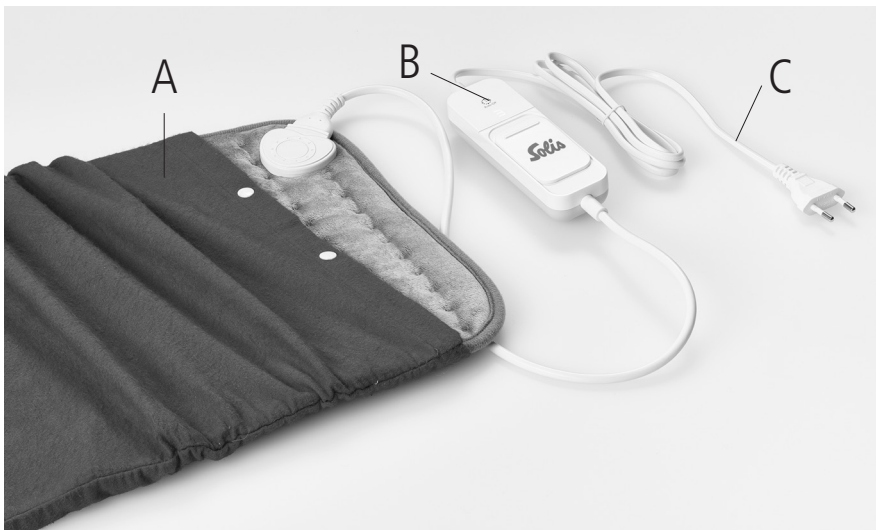Read through this user manual thoroughly before using your appliance

Do not pierce the heating pad, or bring sharp objects near it

Never operate the heating pad when it is folded up or crumpled

This appliance must not be used for children younger than 3

Only use heating pad indoors

Do not chlorinate or bleach the heating pad or bring it into contact with chlorine

Do not iron the pad. The cover can be ironed at a low temperature.

Do not clean with chemicals

Pad and cover not suitable for tumble drying

The heating pad can be either cleaned manually or washed in the washing machine, at a delicate cycle and a max. temperature of 30°C. Please see the label or the washing instructions in this manual

**GB**

## ▲ **IMPORTANT SAFETY PRECAUTIONS**

**GB**

Please read through this user manual thoroughly before using your heating pad, so that you can become familiar with your appliance and can use it safely. We thoroughly recommend that you keep these instructions and, if you pass the appliance on, that you also hand these instructions to the new owner.

For your own safety, please observe the following safety measures when using electrical appliances:

1. To avoid electric shock, neither the heating pad, power cable, control, power adaptor nor plug are allowed to come into contact with water or other fluids. Also, when washing, make sure that the control has been removed along with the power cable!
2. Do not operate the appliance resp. pull or connect the plug with wet or moist hands or feet.

3. Before using the appliance, check that the voltage stated on the type plate of your appliance matches your mains voltage. We do not recommend operating the appliance with a multiway plug adapter.
4. Check your appliance before every use to make sure that the appliance itself, the control, the power cord and the power adaptor/plug are without damage and no parts have come loose.
5. To prevent electric shock, do not use the appliance if the cable, control unit, power plug, the appliance itself or the power adaptor are damaged or if the appliance has other faults, e.g. if the flow of heat is interrupted, the heating pad has fallen down resp. is damaged in any other way or if parts have come loose. Never attempt repairs yourself, but take your heating pad to Solis or a service centre authorized by Solis in order to get it inspected, repaired or mechanically and electrically maintained, **in order to avoid a hazard.**

**GB**

6. Do not leave the power cable hanging over table edges to prevent the appliance from being pulled down. Make sure that neither the appliance nor the power cord or plug comes into contact with hot surfaces such as hotplates or a radiator.

**GB**

7. Never carry the appliance by the power cord and do not place anything on the power cord. Never pull on the power cord to disconnect the appliance from the power mains. Only pull on the plug.

8. Do not fold or clamp the power cord.

9. During use the control is not allowed to lie on or under the pad or be otherwise covered.

10. Never lay the pad on hot surfaces or bring the pad into contact with them.

11. Do not use the appliance outdoors.

12. Unplug the appliance when it is not in use or unattended. Never leave the heating pad switched on unattended. We recommend the use of a

residual current device (RCD safety switch) to ensure additional protection during the use of electrical appliances. It is advisable to use a safety switch with a rated residual operating current of 30 mA maximum. For professional advice consult an electrician.

13. Make no changes to the appliance, control, cable, plug or power adaptor.
14. Only use the heating pad with the control and protective cover provided.
15. Children are not allowed to play with the pad or use it unsupervised. The pad must not be used for small children under the age of three, for children over the age of three, general supervision by an adult is necessary during use, a low temperature setting should be selected. We recommend that children under 8 years of age should generally not operate the heating pad. Pets should not come into contact with the appliance.

**GB**

**GB**

16. People with limited physical, sensorial or mental capabilities as well as unknowing or inexperienced people and children must not use the heating pad, except when supervised and thoroughly instructed on how to use the appliance by a person responsible for their safety. Additionally, they must be aware of the dangers arising from the appliance and know how to operate it safely. Cleaning and maintenance of the appliance must not be done by unsupervised children. Under no circumstances must the appliance be used by children under the age of 3 and be operated by children under the age of 8.
17. Incorrect use may cause injury!
18. Never reach for an appliance that has fallen into water or other liquids. Always use dry rubber gloves to pull the plug before taking it out. Do not start using it again before having it checked for functionality and safety by Solis or a service centre authorized by Solis.

19. Never place the heating pad such that it could fall into water (e.g. next to a bathtub, swimming pool, whirlpool etc.). Danger of electric shock!
20. Never switch the pad on in a damp environment (e.g. bathroom) or when the protective cover is damp. The appliance must always be completely dry before you switch it on.
21. During thunderstorms, disconnect the appliance from the mains and do not use it.
22. Never use the appliance in rooms with flammable or explosive materials or liquids such as, for example, deodorant or hairspray.
23. Never put safety pins or other sharp objects in the heating pad.
24. Never bring sharp or pointed objects or tools near the heating pad. The appliance could become damaged!
25. Never fold or crease the heating pad during use.

**GB**

26. Only operate the appliance as described in this manual. Only use the control provided. Incorrect accessories or misuse can lead to damage of the appliance, fire, electric shock or personal damage.

**GB**

27. This appliance is only for domestic use and is not suitable for commercial use. Not for use at doctors surgeries or hospitals.

28. People with a pacemaker should not use this appliance!

Observe the following when using the heating pad:

- Should you be unsure or need advice on whether the use of a heating pad is good for your health, please consult your doctor.
- Should you notice any continuing pain in your muscles and joints, you should also consult your doctor. These could be signs of serious illness.
- Do not use the heating pad on swollen, inflamed or injured skin areas or on open wounds.

- Should you feel unwell or ill while using the pad, switch off the appliance immediately and set it aside.
- The electronic elements in the control (switch) lead to a slight warming of the control (switch) during use. The control must never lie on or underneath the pad and must never be covered by blankets, pillows or similar items. Nor should you put anything on it, lie down or sit on it.
- Do not use the heating pad on an adjustable bed.
- **Make sure not to sleep on the heating pad.**
- Never switch the appliance off by pulling the plug. Always set the sliding switch to OFF setting at the very bottom first, then pull the plug.
- When using temperature settings 2 and 3, do not use the pad for longer than 2 hours. Prolonged use could lead to burnt skin or cause overheating and damage your appliance. Also,

**GB**

**GB**

you should allow the appliance to cool down fully before switching it on again. If the appliance overheats, switch it off immediately, unplug it and wait until it has cooled down.

- For your safety, the heating pad is fitted with overheating protection. A reliable electronic controller ensures that the temperature does not exceed each chosen heat level and thus prevents overheating of the appliance.

The heating pad should not be used rest. operated by:

- People with mental, physical or psychological restrictions, who cannot use or turn off the appliance independently at all times
- Children under the age of 8
- People under the influence of alcohol or drugs
- People which have sensitive skin and cannot detect warmth or heat
- Sleeping people


**GB**

## ▲ **STARTING UP**

**GB**

- You can use the heating pad to warm up and relax. It is equipped with a sliding temperature switch, allowing you to select one out of three heating settings to find a temperature that feels nice to you. The comfortable and cuddly protective cover must be pulled over the heating element every time it is used and be fastened with the push buttons.
- Unfold the heating pad and make sure that the protective cover is pulled over it and that the connection adapter is firmly connected to the socket in the heating element. Plug the power plug into a suitable power outlet.

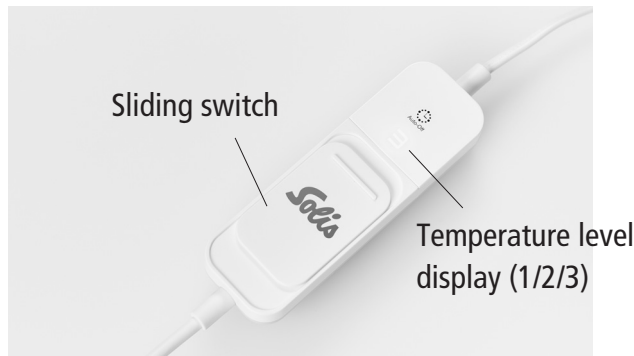- Push the sliding switch up one level, to the position 1. The figure "1" lights up in the display. You can feel the heating pad warm up after a few minutes.
- If you want a higher temperature setting, push the switch up to the next position. The figure "2" lights up in the display. If you want to select the maximum temperature setting, push the switch to the uppermost position so that "3" lights up in the display.
- Should the heating pad get too warm, push the switch back down to the "2" or "1" position. The pad will noticeably cool down after a few minutes.
- If you want to switch the appliance off, push the switch all the way down, the figure in the display goes out. Disconnect the power plug from the power socket.
- **Note:** The heating pad reaches the set temperature only after approx. 15 minutes. Therefore, check the temperature after 15 minutes with the back of your hand to make sure that it feels

**GB**

**GB**

comfortable to you. If necessary, set the sliding switch to a higher or lower position.

- **Caution:** Prolonged use on the highest temperature setting "3" can cause burns on sensitive skin! If you want to use the heating pad for longer periods of time, we recommend using the temperature setting "1".

## **AUTOMATIC SHUTDOWN**

After approx. 90 minutes, the appliance switches off automatically for your safety. If you want to switch the appliance back on afterwards to continue its use, you have to push the sliding switch to the lowest position first, switching the appliance off. The figure disappears. Wait for approx. 5 seconds. After that, you can switch the appliance back on by pushing the sliding switch back to the desired temperature level.

- In any case, the power plug has to be disconnected if you do not want to use the appliance and after you have switched it off.

## ■ OVERHEATING PROTECTION

For your safety, the heating pad is equipped with overheating protection. An electronic control device makes sure that the temperature of the selected temperature level is not exceeded, preventing the appliance from overheating.

**GB**

## ▲ STORAGE AND CARE

- Always let the heating pad cool down completely. **Caution:** The heating pad has to be completely dry and cool before you fold and/or store it.
- Do not place anything on the heating pad while it is in storage. The storage space should be dry, clean and not exposed to direct sunlight.

## ▲ **WASHING**

- The cover and the heating element can be washed on the delicate cycle of the washing machine at 30°C or by hand. Wash the cover separately or with similar colours as it might stain. To do that, please proceed as follows:

**GB**

1. Before washing, the power plug has to be disconnected from the power socket. Let the heating pad cool down for at least 10 minutes.
2. Now, disconnect the connection adapter from the socket on the heating element and take the power cord with the control element off the heating pad.

**Caution:** Power cord, control element, connection adapter and power cord must never come into contact with water!

3. Wash the cover and the heating element either on the delicate cycle of the washing machine at a max. temperature of 30°C or hand-wash both items in lukewarm water. To do that, place the heating element and the cover in the bathtub, add lukewarm water and a small amount of mild detergent and carefully squeeze it out. **Caution:** Never crease or wring out! After that, rinse both items repeatedly in lukewarm water until all detergent residues are washed out. Leave both parts to air dry, lying flat on a towel. No items must be placed on the heating element or the cover.
4. When the heating element and the cover are completely dry, pull the cover back over the heating element and connect the connection adapter to the socket. You can now start using the appliance again or put it into storage.

**GB**

## ▲ TECHNICAL INFORMATION

Model-No.	HP18303, Type 2260
Product name	Heating Pad
Voltage/Frequency	220-240V~ 50 Hz
<b>GB</b> Power	100 Watt
Dimensions	40 x 30 cm
Weight	0,375 kg

Technical changes reserved.

## ▲ DISPOSAL

EU 2012/19/EU

Information for correct disposal of the product in accordance with the European Directive EU 2012/19/EU.

**GB**

At the end of its working life this equipment must not be disposed of as an household waste. It must be taken to special local community waste collection centres or to a dealer providing this service. Disposing of electrical and electronic equipment separately avoids possible negative effects on the environment and human health deriving from an inappropriate disposal and enables its components to be recovered and recycled to obtain significant savings in energy and resources. In order to underline the duty to dispose of this equipment separately, the product is marked with a crossed-out dustbin.

Solis of Switzerland Ltd reserves the right to carry out technical and visual changes and modifications to improve the product at any time.

**GB**

## **SOLIS HELPLINE**

This item stands out due to its long-life cycle and reliability. Despite this though, if a malfunction ever occurs, simply give us a call. Often a breakdown can be easily and quickly remedied with the right tip or technique without having to straight away send off the device for repair. We will give you all the advice and practical support you require. The contact addresses can be found at the end of this manual.

Keep the original packaging in a safe place so that you may use it in the event that transportation or shipping of the product becomes necessary.

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## 2 YEAR GUARANTEE

(Only valid with a till receipt)

We provide a guarantee for this product, if the defect has demonstrably arisen as a result of faults in the material or construction, and has arisen despite proper handling and care. The guarantee starts from the date of sale, as noted on the sales receipt, which must be enclosed with the item. For commercial use, the guarantee lasts for 12 months. The detailed guarantee conditions are available at [www.solis.com](http://www.solis.com).

### CUSTOMER SERVICES

It is always worth calling us before sending in the product, as malfunctions can often be easily solved by the appropriate tip or trick from our experts. Please visit us at [www.solis.com](http://www.solis.com) for country-specific contact details. Naturally, we would also be pleased to provide you with information over the phone. Please contact us via the following telephone numbers:

For Switzerland: +41 91 802 90 10  
For Germany: 0800 724 0702  
For Austria: 0800 22 03 92  
For the Netherlands: +31 85 4010 722  
For Belgium: +32 2 808 35 82  
For the other countries: +41 44 874 64 84

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